



WOXRf-P

Microneedle RF Machine




WOXRF-P

Microneedle RF Machine


The WOXRF-P is an advanced microneedle RF device designed for versatile skin treatments. Equipped with dual treatment modes and an adjustable needle depth system, it ensures precision and effectiveness. The product offers two handpiece options: the Microneedle Handpiece and the Dual Max SRF Handpiece, catering to facial and body applications. It employs a sophisticated needle shooting system with interchangeable needle matrices (MRF-49 and MRF-25) for tailored results. Ideal for skin rejuvenation and tightening, the WOXRF-P enhances treatment flexibility and outcomes.

How Does It Work? Microneedle Handpiece

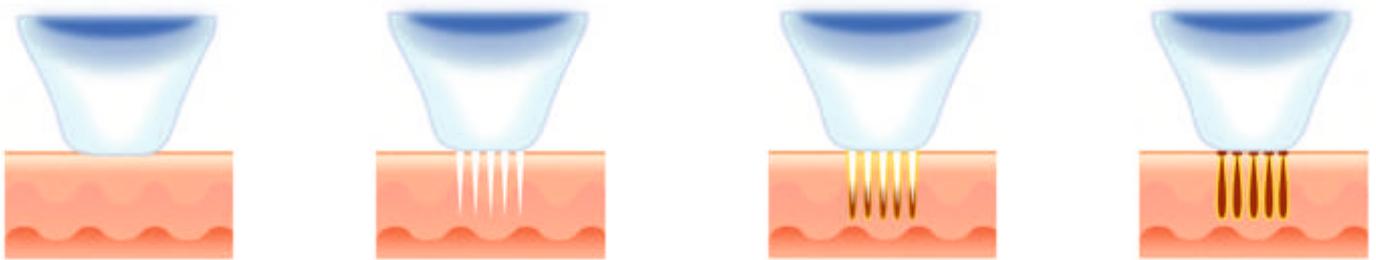

Principle of Treatment : Using a professional radio frequency micro needling device with fine needles to create small pores in the surface layer of the skin. It creates a triggering effect that promotes the natural production of new elastin and collagen. Radiofrequency delivers RF energy deep into the dermis to increase the impact of traditional micro needling.

Dual Max SRF Handpiece

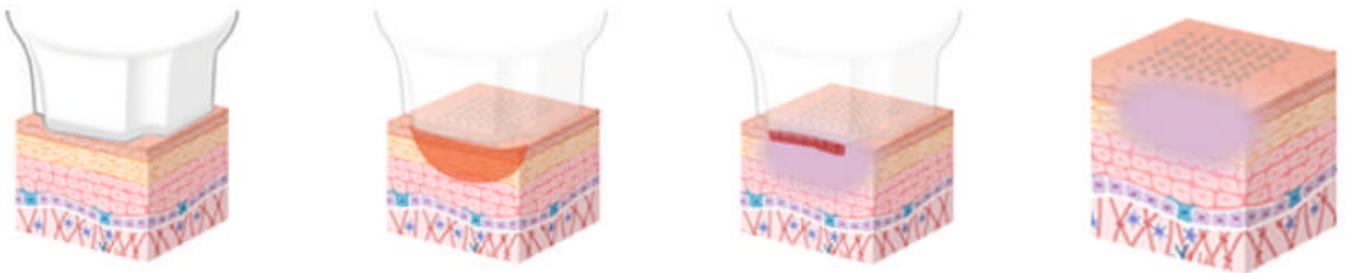

Principle of Treatment : The patented flow needling technology serially inserts the microneedles row-by-row for a more comfortable procedure. Serial insertion reduces unnatural stretching in order to further reduce pain and downtime. The depth of the coagulation zone is precisely controlled. Sub-dermal delivery of RF energy eliminates thermal damage to the epidermis, providing faster healing time for your patients. The adjustable depth of the RF energy delivery allows for multiple passes and treatment of delicate areas.

Key Advantages

- **Two Treatment Methods:** There are two treatment methods with a dual matrix tip and RF microneedle tip to meet the needs of different clients.
- **Needle Depth Adjustable :** The needle is 0.3~3 mm adjustable, operates epidermis layer and dermis layer by controlling the needle depth in units of 0.1 mm.
- **Needles Shooting System :** Automatic output controlling, can make the RF energy better distributed in dermis, so patients can get better treatment results.

Technical Specifications

RF Frequency	5 MHZ
RF Energy Intensity	1~10 Level
Power	80W
Type of Needles	81 tips, 49 tips, 25 tips
Depth of Needle	0.3-3mm(Adjustable)
MRF Head Area(cm2)	1*1, 1.5*1.5, 2*2
SRF Head Area	64 pin / 2*2 cm2
Input Voltage	110/220V; 50/60HZ

Microneedle Handpiece

Sterilized Disposable Needle Cartridges



3 Size Cartridges Suits To Different Treatment Area



MRF-81 Needles



MRF-49 Needles



MRF-25 Needles

Dual Max SRF Handpiece

Dual Max SRF Handpiece



SRF




Treatment Results




Before




After

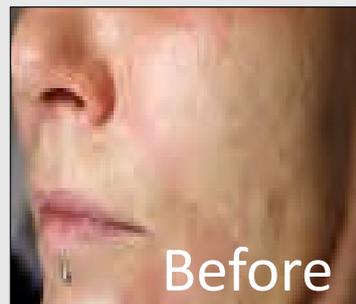

Before

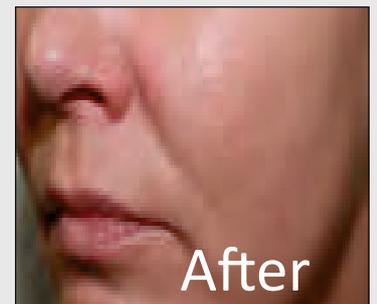

After

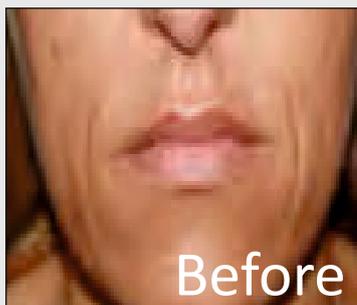

Before



After



Before



After